

Sendoline®

PERFECT ENDO



S5

ROTARY SYSTEM



S5 ROTARY

S5 Endo Motor

SAFE, SIMPLE AND EFFICIENT!

Sendoline's rotary system for mechanical cleaning.

Features:

- 1. Easy to use**
S5 Endo Motor + S5 Rotary Files STERILE = S5 Rotary System.
Uses only 5 settings with 5 files (S5 Rotary Files STERILE)
Complete versatility, works with most rotary file systems
Cordless freedom: no cables
- 2. Security**
2 options: auto reverse / auto stop
Torque indicator: Light indicates when 75% of the torque is reached
Low speed: 300 rpm
- 3. Reduced treatment times**
Treatment can be completed safely, and more efficiently using the S5 Rotary System with the unique S-Profile
- 4. Accessibility**
Only 11 mm diameter contra angle head
Shorter shanks (13mm) on rotary files
- 5. Non-stop**
2,5 hrs continuous use
- 6. Ergonomics**
Weight: 120g
- 7. Apex Clip**
Easy connection to an apex locator


Apex Clip

Apex Clip for easy connection to an apex locator.

Safe and secure

Light indicates when 75% of the torque is reached

5 settings = 5 files

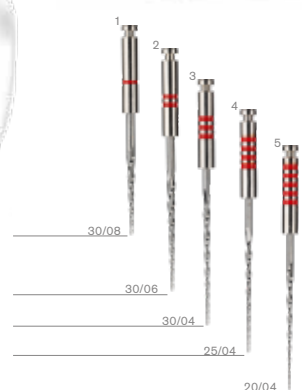

www.sendoline.com

SYSTEM

S5 Rotary Files STERILE



1, 2, 3, 4, 5 – FIVE STEPS TO PERFECT ENDO.

The S5 Rotary System offers a simple and effective solution for the dentist.

Features:

1. Only 5 files and all with the unique S-cross section
2. Shorter shank (13mm) for better access
3. Long progressive flutes to prevent the “screw-in-effect”
4. Enhanced debris transportation
5. Presented sterile in a blister pack

30/08 30/06 30/04 24/04 20/04
ISO/TAPER

“The S5 file system is so simple to use, it gives me predictable results every time!” Dr. Paul Godbold BDS General practitioner, UK

*“The S5 rotary system helps me to efficiently and safely create a canal shape that is perfect for cleaning, disinfecting and then obturating every canal I come across.”
Dr. Claire Redley BDS General Practitioner, UK*

S5 Apex Locator

A highly accurate electronic apex locator which allows precise measurement of the root canal’s working length. Its compact design makes it portable and easy to use.

Features:

1. High precision of apical foramen localization in wet and dry canals
2. Fully automatic measurements
3. Compact and user-friendly design
4. Automatic turn-off
5. Audio feedback with volume control



Preparation

Big Apical Files STERILE

For the best results in cases of teeth with larger apical diameters use Big Apical Files STERILE, a perfect completion for the S5 Rotary System. The markings on the Big Apical Files STERILE correspond to the program (3) of the S5 Endo Motor. Shape the canal as usual with the S5 Rotary System and simply finish the apical portion using the Big Apical Files STERILE. Available in sizes: ISO/TAPER 35/04 - 50/04 and file length 23 or 28 mm.



Obturation

S5 Gutta Percha Points

The S5 Gutta Percha Points are specifically designed to match the canal shape created with the S5 Rotary Files STERILE. Finish the preparation with either file no. 2 or no. 3 and simply select the S5 Gutta Percha Point corresponding to that file. Package: 60 pcs/box.



S5 Paper Points

The S5 Paper Points are pre-sterilized and specifically designed to match the S5 Rotary Files STERILE. They are hand rolled to exact dimensions to match the canal shape created with S5 Rotary Files STERILE size 2 and 3. Available in sizes: ISO/TAPER: 2 (30/06) och 3 (30/04). Package: 3 pcs/blister, 120 pcs/box.



Retreatment

S5 Retreatment Files

The S5 Retreatment Files enables the removal of old filling materials quickly and easily. To conduct the retreatment, simply apply the Crown Down technique using only 4 files. The markings on the S5 Retreatment Files correspond to the programs (1-4) of the S5 Endo Motor. ISO/TAPER file 1: 12/30 (17 mm) and file 2-4: 08/30, 06/30 and 04/30 (23 mm).



Sendoline[®]
PERFECT ENDO

Sendoline AB
Box 7036, SE-187 11 Täby | Måttbandsvägen 12B, SE-187 66 Täby | Sweden
contact@sendoline.com | www.sendoline.com

The contents of this brochure including layout and images are protected by copyright. Sendoline, S1, S5, S5 Rotary Files, S-Files / S-fil, including the S symbol and NiTi-S, are Sendoline brands and/or images. ©2020 Sendoline AB. All rights reserved.